

A woman with long brown hair, wearing a white sleeveless dress, is shown in profile, looking upwards and to the right. She is reaching out with her right hand towards a large, glowing white circular interface that appears to be part of a futuristic machine. The background is a light, textured wall with vertical lines.

SYNETICA

EVOLUTION

The future is now



The maximum evolution
of body contouring systems

A completely Italian
Design



CHARACTERISTICS:

- Weight: 27 Kg
- Methodic: ultrasound cavitation with radiofrequency
- Max output power cavitation: 40 Watt
- Cavitation frequency: 40 kHz
- Max output power radiofrequency: 25 Watt
- Frequency radiofrequency: 470 kHz
- Handpieces: 2 for external treatments, 3 probes for internal treatments
- Software upgrading: USB, WiFi



SYNETICA Evolution: Internal Lipo Cavitation

An extraordinary solution for body contouring has been created that promises to be decisive for aesthetic problems related to fat deposits: Synetica Evolution. Synetica Evolution is the lipo cavitation of the future: compact, effective and unique – because of its characteristics – among the methods used by cosmetic doctors.

The Evolution of Cavitation


Synetica Evolution is the most advanced development of the cavitation concept, because it acts from within and targets the most restricted areas that, so far, are untreatable by the cavitation equipment on the market. Its enormous potential is concentrated in the "needle" hand-piece, which transmits the contouring action directly inside the tissue, aiming it at the most affected points.

Internal Action

With minimum invasiveness, the cannula is inserted under the skin area to be treated. The direct solicitation of fat greatly enhances the effect of dissolution of fat masses, because the cavitation action stimulates fatty tissue in a more immediate way. However, the effectiveness of Synetica Evolution is even greater, because the ultrasonic waves are generated on the entire length of the inserted cannula, not only at the tip.

Micro Action


Synetica Evolution cannulas have a thin diameter of 2 to 3 mm and penetrate areas that are normally the most difficult to reach. For this reason, it is the method used to counteract the fatty deposits in areas hitherto impossible to treat, such as chin, armpit, inner thigh and knee and ankle areas. Moreover, thanks to the precision of the cannula's action, it is also possible to work in close proximity to the most delicate areas, without the risk of damaging them.

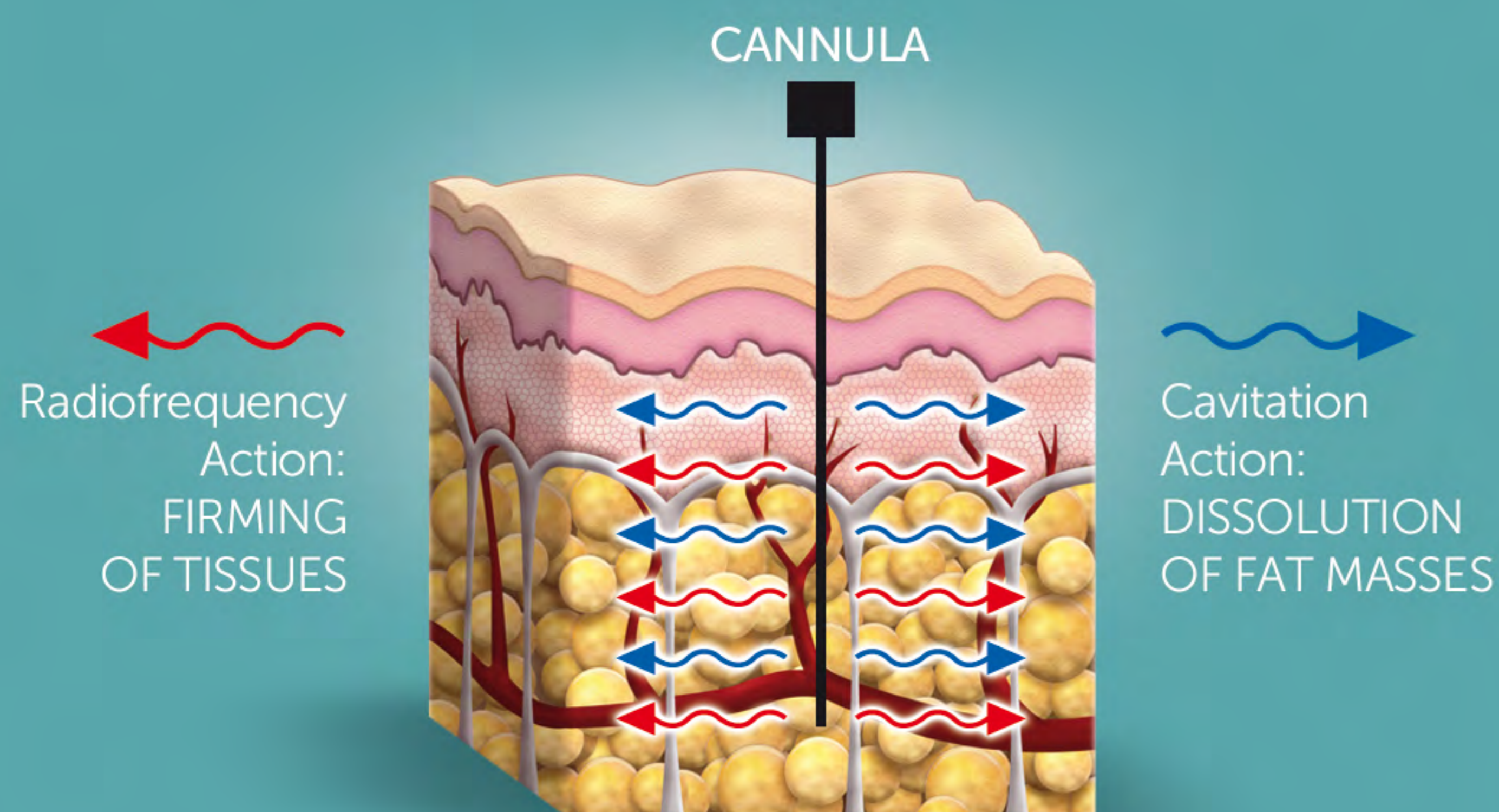

SURGICAL PROCEDURE


1

Percutaneous anesthesia

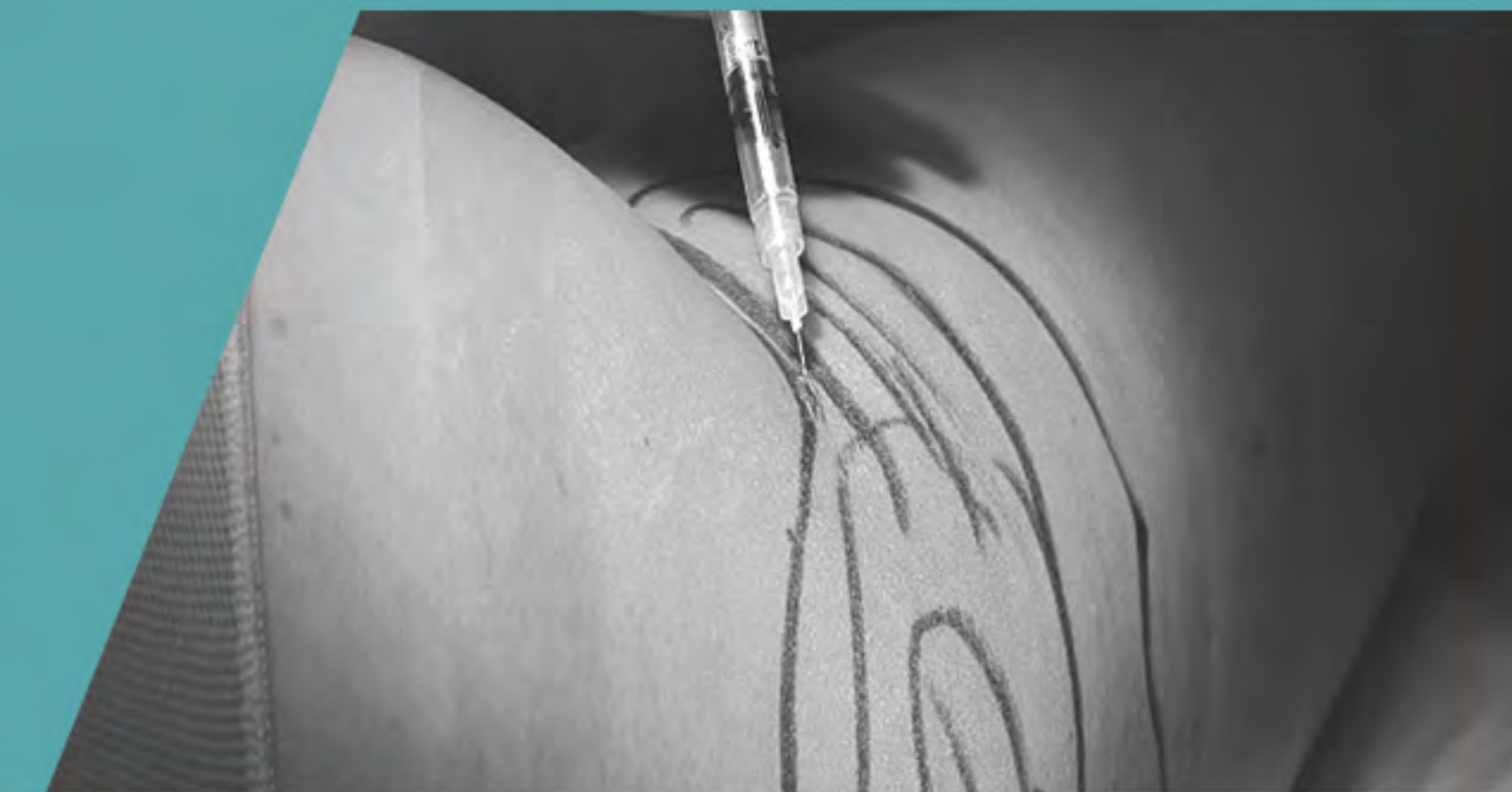

2

Tumescent infiltration

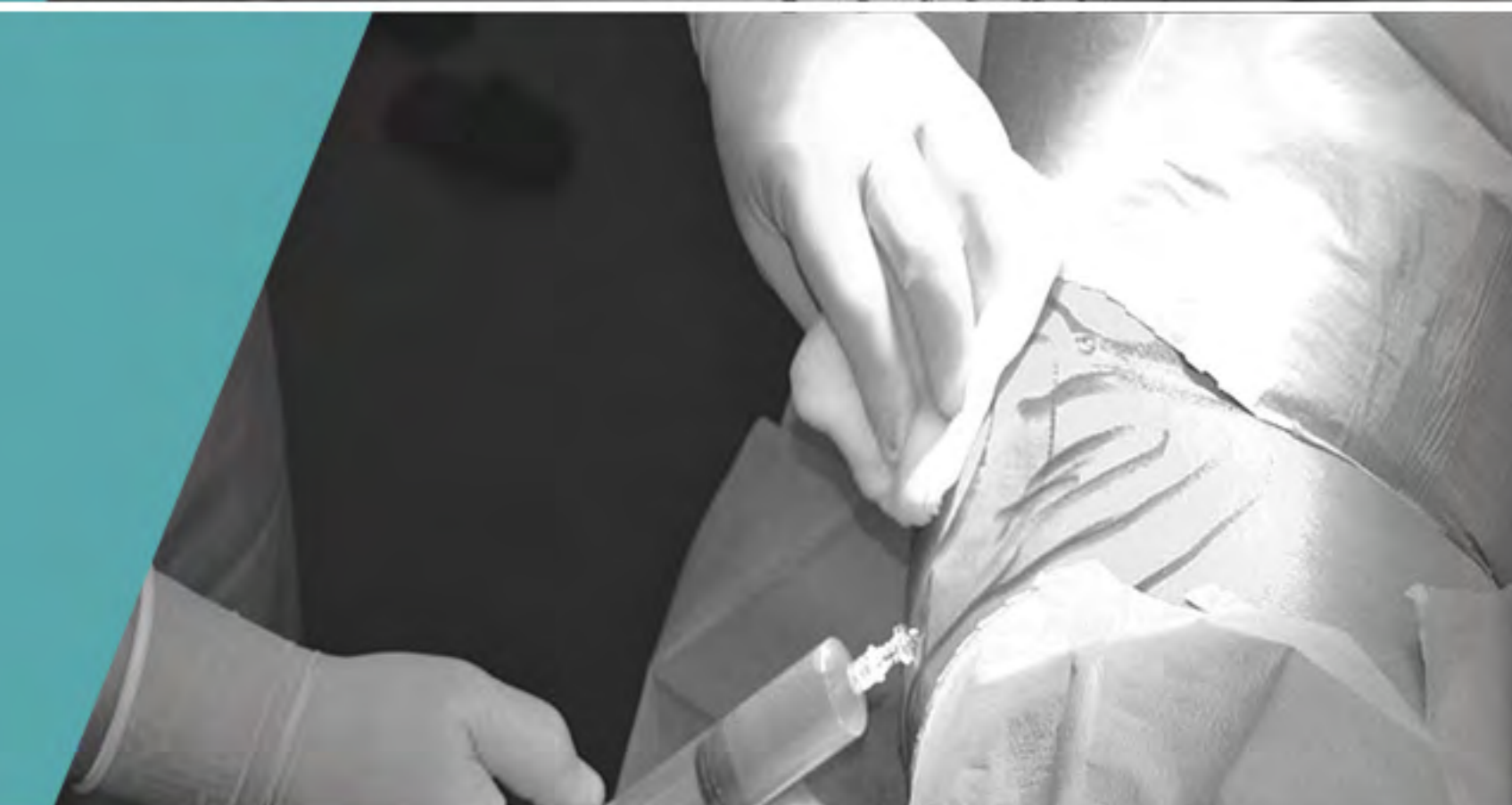


3

Lipo cavitation



4

Squeezing

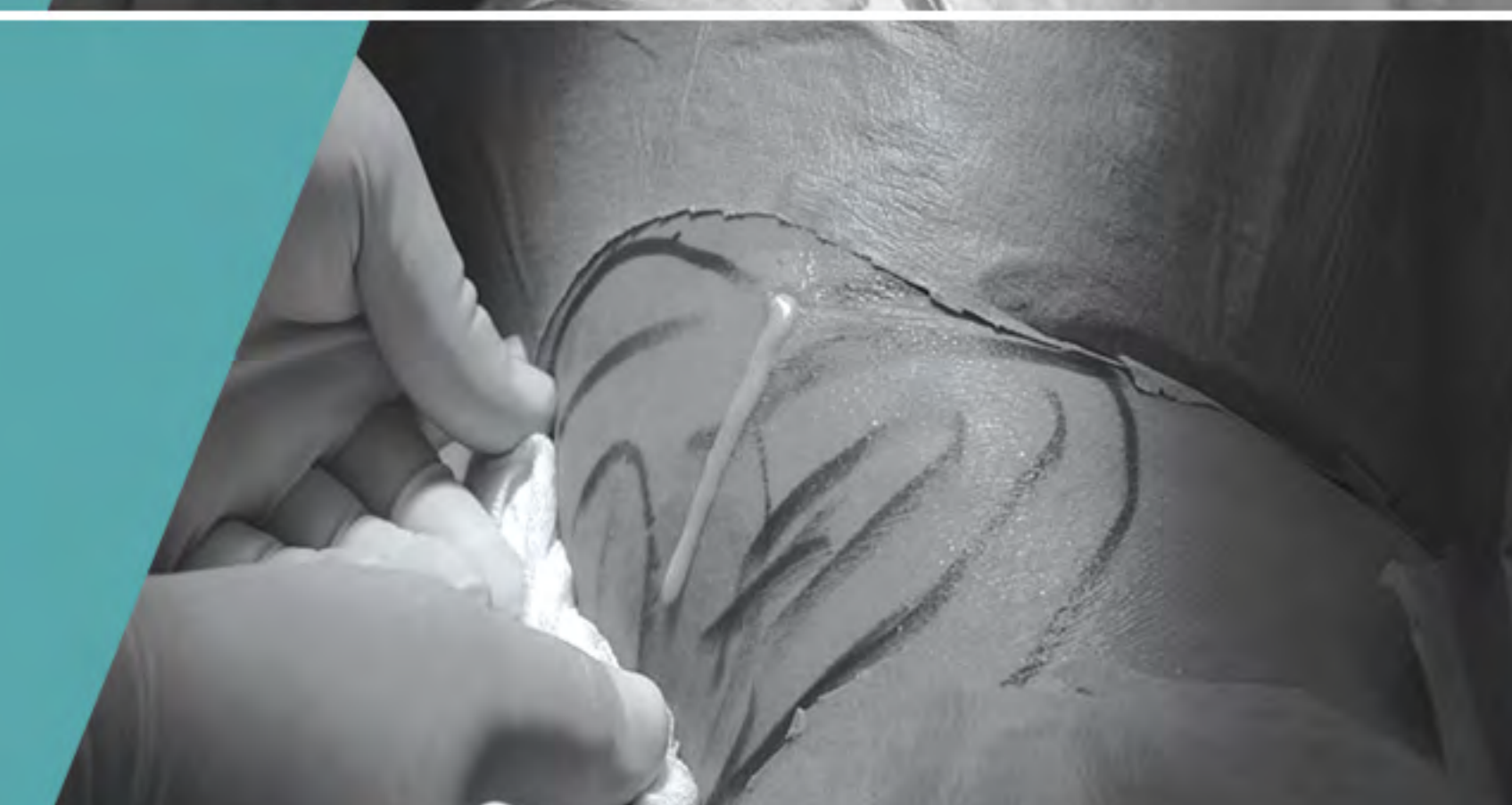

5

Lipo Suction



6

Fats collection



Applications

LIPOSCULPTURE TREATMENTS IDEAL FOR:

- Chin
- Arms
- Knees
- Ankles
- Hips
- Buttocks
- Abdomen
- Inner and outer thighs

INTERNAL ACTION with lipo cavitation and radio frequency

The great modernity of Synetica Evolution lies in uniting the potential of internal lipo cavitation and the firming effect of internal radio frequency. In this way, there is no risk of replacing unsightly fat mass with the flaccidity of treated tissues, because action is taken on both fronts simultaneously.

Radio frequency waves are simultaneously emitted from the body of the cannula, compensating the dissolution and depletion of tissues through micro cavitation. Radio frequency is also delivered along the entire length of the cannula and not only in the terminal part. In this way, the effect of skin rejuvenation induced by radio frequency is continuous, widespread and, above all, parallel to cavitational action.

Through the thin body of the cannula, radio frequency reaches the smallest gaps treated, in such a way as to leave no sign of depletion due to cavitation action. As the cannula passes, firming action takes place to leave tissues firm immediately and not oblige the patient to resort to subsequent contouring treatments. In this way, the patient remains completely satisfied.

Results





Treatment with Synetica Evolution

The biggest advantage of equipment such as Synetica Evolution lies in offering liposuction surgery results without subjecting the patient to the stress of hospitalisation or general anaesthesia which is often poorly tolerated. In fact, treatments take place in an outpatient setting and, when they are over, patients are fully able to resume their routines.

The treatment, performed by a medical specialist, does not require an operating room or sedation. Only the part to be treated, where the cannula to be used will be inserted, is anaesthetised. At the end of the session, patients can resume their daily activities because there are practically no contraindications: in fact, the incisions made for inserting the cannula are small (2 mm) and are closed with a simple butterfly bandage.

The results obtained are already visible at the end of the first treatment; in fact, it takes just one session to eliminate the fatty deposits. This is a real liposculpture thanks to its noticeable effects, but of the soft type due to the minimum of invasiveness required. The dissolved fatty tissue is aspirated on the spot and, in part, disposed of by the body itself.

A safe choice

Synetica Evolution is an extremely avant-garde system for the dissolution of fatty deposits and PEFS. Approved by cosmetic doctors, it is equipment that is particularly appreciated because it works in a soft way to obtain significant results. The response of patients is also excellent, because they can return to normal daily activities immediately.

Advantages of Synetica Evolution

EXCLUSIVITY: Synetica Evolution is the only equipment with well-developed characteristics on the market for body contouring aesthetic equipment.

IMMEDIACY: the results obtained with Synetica Evolution treatment are immediately visible; the first treatment already provides tangible proof of the effectiveness of the equipment.

ACCURACY: thanks to the possibility of operating with thin cannulas, Synetica Evolution permits targeted treatment, even in extremely restricted areas that were previously untreatable.

OUTPATIENT TREATMENT: treatment with Synetica Evolution does not require an operating theatre or anaesthesia, but can be carried out as an outpatient.

SATISFACTION: Synetica Evolution fully satisfies patients' demands, because it reconciles depletion of fatty deposits with firming of emptied tissues.



Al Shirawi
Healthcare Solutions
Solutions You Can Trust

Office: Building No. 26, St - 8A,
Industrial Area 2, Ras Al Khor,
PO Box 93, Dubai, UAE
Contact: +971-4-3337644
Email: info@alshirawi.health
Website: www.alshirawi.health



The Intelligent Beauty



synetica-evo.com

3W SRL

Via Don Giuseppe Mucciardi 5 (Zona Industriale)
86020 Campochiaro (Cb) - ITALY

Tel.: +39 0874 1896435

Fax: +39 0874 1892046

Email: info@triworks.it

www.triworks.it



ASSOBIO MEDICA

Is the Italian National Federation which represents companies operating in the field of biomedical technology, diagnostics, medical devices, and telemedicine